BEAUTIFLY

B-Derma Ice
PRO
series

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Technical data:

Product size:	190*45*40mm
Product weight:	155g
Working voltage:	3.7V
Charging voltage:	5V
Material:	PC+ABS
Suction power:	61Kpa±5Kpa
Battery capacity:	1000mAh
Charging Time:	150min

Part description - device

1. Replaceable head
2. On / off button intensity control
3. Hot and cold compress/switch key
Hot and cold compress
4. Charging port
5. Air inlet

Part description - heads**Microcrystalline tip**

It can be used on the whole face with the exception of the eyes and in areas on which hard skin such as the knees or elbows were formed.

**Big round hole tip**


Except eyes. Can be used on any part of the face.

**Small round hole tip**

Can be used on sensitive skin around the eyes, and also other parts of the face.

**Oval hole tip**

Can be used at eye corner, nose-wing and skins around the mouth.

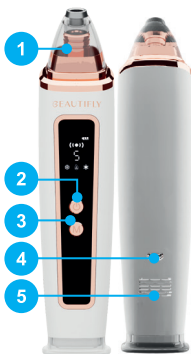Principles of operation

Lifting massage

The vacuum massage system improves the firmness and oval of the face. Vacuum massage supports better blood and lymph circulation, making the skin tense, better supplied with blood, and reduced swelling and sagging skin. The technology of vacuum absorption, used for skin lifting and massage, increases the absorption of cosmetics. It makes the complexion brighter and radiant. At the same time, it strengthens the skin's defenses and keeps it supple and elastic.

Microcrystalline head

Produces over 100,000 natural abrasive particles that remove dead skin cells from the face. The complexion remains smooth and delicate.

Blue photon light

Intended for acne-prone, but also sensitive skin. Bio-effects: has a positive effect on acne skin, reduces inflammation, regulates the work of sebaceous glands and the amount of sebum secreted, the skin becomes tighter and enlarged pores become less visible.

Red photon light

It works at the level of the epidermis and dermis. It reduces and soothes all kinds of irritation. In addition, red light improves skin elasticity and elasticity. It delays the aging process and reduces wrinkles and discoloration. It is also great at reducing stretch marks, scars and cellulite.

Green photon light

It works at the level of the epidermis. Recommended for discoloration and sun spots. It penetrates into the basic layers of the epidermis, reduces expression lines, stimulates nutrients, and helps regenerate mature skin.

Area of use

B-Derma ICE is equipped with 4 types of heads, each of which has a different application.

Microcrystalline head: can be used on any part of the face except near the eyes. It can also be applied to areas of the body covered with hard calloused skin, e.g. on the elbows or knees.

Before using the massage heads, it is necessary to use a cosmetic that gives a glide on the skin, e.g. cream, oil or serum.

Large round head: can be used on any part of the face, in addition to the area around the eyes. Features: Strong suction, suitable for larger surfaces (e.g. forehead). While using, slowly move over the skin and do not stay in one place. Pan from bottom to top as you move slowly from chin up.

Small round head: can be used on sensitive skin around the eyes, also on other parts of the face. Features: Weak suction power. It can be used for areas where the skin is thin and delicate. Not suitable for heavy massage.

Oval head: can be used on any part of the face and neck, except close to the eyes. Features: It is intended for lifting and firming, improves blood circulation in the skin and increases its elasticity.

Usage

1. On / Off

Press and hold the power button  for 3 seconds, the product is on, the display is on, and the suction power shows level 1. Then press the same button to start the next suction levels.

2. Power adjustment

Briefly press the switch button  to adjust the suction power, each time you press it to increase one level and so on (1-2-3-4-5-0-1-2-3-4- 5)

3. HOT mode

Briefly press the  button under the on state. The  HOT lamp will turn on and automatically switch to HOT mode. Place the HOT surface on your face and move it over the entire surface in circular motions to loosen the skin and open pores.

4. ICE mode

Press the  key twice under the on state. The lamp  in ICE mode will turn on and switch to ICE mode automatically. Place the ICE surface on your face and move it around in circular movements to soothe the skin, close the pores and regenerate it after the microdermabrasion treatment.

The device heats up after using the ICE mode for a long time, so it is not recommended to use this mode for more than 3-5 minutes. After use, turn off the product and wait 30 minutes before using it again.

5. Photon light therapy

Press the **M** key three times under the on state to enter the photon light mode. With this feature enabled, move the device across your face in circular motions. Inside out, bottom up, symmetrically left and right. Treat each area for about 30 seconds.

6. Precautions

Set the suction head to 90 degrees and slowly move the suction head over your face. Do not hold in the same part for a long time. (Before removing the suction head from the skin, tilt it to the side to release the air under pressure. Then you can pull it off safely without irritating the skin. Do not pull or move it directly to avoid pulling the skin and redness).

Application

- Discoloration
- Scars
- Stretch marks
- Extended, very visible pores
- Comedones and acne (avoid inflammatory changes, phlegias)
- Wrinkles
- Skin devoid of a healthy color with the first signs of photoaging

How to use

Microdermabrasion

1. Please thoroughly wash your face. This will makes your pores open up.
2. Please do not apply the device at the same spot on your skin more than 3 seconds. It needs to move up and down constantly. We recommend that you include photon blue light, which is antibacterial.
3. Please attention that sunscree moisturizing in time to avoid skin hurt after using the product, it is recommended to apply toner or cold water on the face, this will close the pores.
4. Please adjust the suction pressure based on your skin type.
5. After used, it is normal that your face to be a little red and sensitive, the redness will disappear in 30 min.
6. After the treatment is over, wash your face with a soothing lotion and use a moisturizing pre-treatment with a filter. We recommend using a moisturizing mask and using the ICE function along with photon light therapy.

The lifting massage

1. Use a cleansing liquid or peel-off mask, and use the HOT mode for about 3-5 minutes to open the pores and loosen the skin.
2. Apply the appropriate massage tip to the device.
3. Turn on the photon light according to the effects you want to achieve.
4. Point the tip at the skin at a 90-degree angle. You can use any suction level that suits your skin's needs.
5. During the massage, you should move the device slowly and each time in the same direction - do not move the device oval. You can suck in and let go of the skin where you need it. Do not use in places of purulent acne or other skin inflammations.
6. Do not use the device on the same parts of the face for too long, redness may appear on the skin. In that case, please stop using for a while
7. After the treatment, ICE mode can be used to soften the skin and close the pores.

Tips:

Do not wash the device with water. The B-Derma ICE device is not waterproof.

Contraindications

- Active acne phlegosis
- Rosacea
- Therapy with retinoids
- Acidic, lasers - at least 4 weeks before microdermabrasion
- Psoriasis
- Virotic, bacterial and fungal diseases of the skin vascular skin
- Vascular skin
- Scarring
- Distinct continuity of the skin (fresh wounds)
- Skin cancer or suspected lesions in the place undergoing surgery
- Fillers in the face - at least 4 weeks before microdermabrasion, otherwise, the preparation may move.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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English

Hereby, Beautifly Sp. z o. o. declares that the device equipment Beautifly B-Derma Ice is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
 2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device.
 3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
 4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
 5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.
- In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@beautifly.eu

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I represent that I have read and accepted the terms and conditions specified in the Warranty Certificate.