BEAUTIFLY

B-Glossy

Blush

series

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B–Glossy specifications

Rated voltage:	5V/1A
Rated power:	4W
Battery capacity:	600mAh
Charging time:	3 hours
Massaging head diameter:	35 mm

Description of the appliance

1. Ultrasound head
2. Photon lights
3. Battery level
4. Light therapy: Red/Blue
5. Operating mode
6. Buttons on/off
7. Button mode


Ultrasound therapy

The frequency of the ultrasound is 3 MHz. The ultrasound wave has a strong power and high energy. Ultrasound has a thermal, physicochemical and mechanical effect. This means an increase of the penetrability of cell membranes and stimulating cell metabolism. During therapy, there occurs a phenomenon in which active substances from cosmetic products are surrounded by transport vortices and forced into the deeper layers of the skin. Deep tissue massage improves circulation, evens skin tone and makes the skin feel smooth.

1.Mechanical operation: based on the use of ultrasound as a micro-massage. Ultrasounds improves blood supply, thus changing the structure of cell membranes and proteins. Due to this, the permeability of skin tissues is increased, the production of collagen and elastin is intensified, which perfectly translates into regeneration processes and regulation of the work of the sebaceous glands.

2.Thermal effect: the ultrasound warms up the skin surface. It improves blood circulation and increases the amount of nutrients in skin cells. It doesn't cause a heat sensation because the heat is produced internally.

3. Chemical operation: ultrasound improves the ability of catalysis, that is an increase in the speed of skin regeneration processes. Better permeability of cell membranes and improved metabolism help to significantly reduce swelling, redness and lividness on the face. With the use of suitable active substances, the visibility of capillaries and erythema become much smaller.

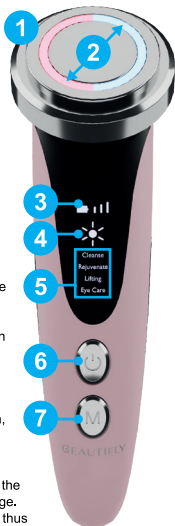Photon light therapy

Red light

It is used for wrinkle removal. The wavelength is $625 \text{ nm} \pm 3 \text{ nm}$ and can be applied to all skin types. Bio-effects: stimulation of the cells to produce collagen and elastin faster, increasing of cell viability, acceleration of blood circulation. It has a firming effect quickly smoothing the skin. It perfectly reduces discoloration evening out the skin tone. It helps remove puffiness and dark circles under the eyes.

Blue light

Blue light with wavelength of $465 \text{ nm} \pm 3 \text{ nm}$ is intended for acne-prone and sensitive skin. Bio-effects: it has a beneficial effect on acne-prone skin, reduces inflammation, regulates the work of sebaceous glands and the amount of sebum secreted. It tightens the skin, making enlarged pores less visible.

Ion+ Positive ion therapy

Consists of cleansing the skin from dead corneum and any other impurities. This improves the permeability to active substances. Discolorations become less visible with time, and the skin is visibly brightened. This effect will be much better if you combine the therapy with, for example, micellar water.

Ion– Negative ion therapy

Consists of the use of a current of sufficient intensity to overcome the barrier that prevents deep absorption of active substances from cosmetic products. After a series of treatments, the skin becomes visibly taut, smoothed, hydrated and brightened – it acquires a radiant glow. This effect will be much better if you use a proper cosmetic with the Ion– therapy, for example serum or ampoule.

Vibration massage

Pleasant massage has a relaxing effect and reduces muscle tension. Stimulating the work of muscles significantly improves the shape of the face.

4 modes of operation – methods of use

Caution:

Begin the treatment with cleansing (Cleanse mode), use a facial cleanser. Thorough cleansing of the pores will involve a deep removal of impurities, which will come out to the surface of the skin. Wash the surface of the skin with makeup remover. Set each mode for 4 minutes of use time. Within 1 minute of inactivity, the appliance will automatically switch off. The appliance works on the basis of induction, which means it starts working when it comes into contact with the skin.

Use of Cleanse mode (with a cotton pad)

In this mode you have the choice of 3 levels of intensity, heat and vibration. Thanks to this, you can effectively cleanse the pores of impurities, brighten the skin and improve its tone.

1. Remove the cotton pad clip.
2. Attach a cotton pad.
3. Pour makeup remover onto the cotton pad.
4. Press and hold the button for 3 seconds to switch on the appliance.
5. Put your fingers around the sensor strips.
6. Cleanse your face with circular motions.

TREATMENT TIME: 4 minutes

Use of Rejuvenate mode (no need to use a cotton pad)

This mode uses red light, heat and vibration massage, as well as blue light, heat and vibration massage. It has a strong firming effect stimulating the skin to intensive production of collagen and elastin. It effectively levels wrinkles and furrows, moisturizes leaving the effect of rejuvenated skin.

1. Apply serum or mask on your skin
2. Press the power button for 3 seconds to switch on the appliance and select Rejuvenate mode.
3. Holding the appliance, touch the sensor strips located on both sides of the appliance. Perform the treatment moving the appliance from the bottom of your face upwards.

TREATMENT TIME: 4 minutes

Use of Lifting mode (no need to use a cotton pad)

This mode effectively stimulates and relaxes the skin, as well as it helps to deeply moisturize it. Red light combined with vibrations and heat is great for remodelling of the face, letting it regain its beautiful youthful oval.

1. Press and hold the power button for 3 seconds to switch on the appliance. Select the Lifting mode.
2. Holding the appliance, touch the sensor strips located on both sides of the appliance. Move the appliance from the bottom of your face upwards.

TREATMENT TIME: 4 minutes

Use of Eye care mode

This mode uses vibrations, heat and optionally red light. After the treatment on your facial skin, we move on to the eye area and do a precise treatment of the delicate skin around the eyes. The skin of the eye area requires special care, and the right treatments will help you level laughter lines, puffiness and dark circles around the eyes.

TREATMENT TIME: 4 minutes

Operating instructions

Use of charger

1. When the symbol  appears on the controller display, it indicates a low battery level. Insert the charger plug into the controller DC jack, next insert the controller plug into a socket, the battery level indicator on the LCD display will start to flash.
2. The rechargeable battery life-span is equal to or greater than 500 charge-discharge cycles in normal conditions, which allows it to be used for over 3 years (once a day).

Precautions

1. Do not leave the appliance plugged in for a long time.
2. Do not touch it with wet hands, it may lead to an electric shock or short circuit.
3. Do not put the appliance in places of high humidity, as this may result in electric shock, short circuit or damage to the insulation.

4. Keep the appliance out of the reach of children.
5. The device must not be used by children.
6. The device must not be used in places where flammable materials are stored.
7. Do not clean the appliance with corrosive substances.
8. Use only the charger included.
9. Make sure the appliance is dry before connecting to the charger.
10. Do not use the appliance on sunburned skin.
11. Stop using in case of undesirable side effects.

Contraindications

- pregnancy
- metal implants, cardiac pacemakers
- metal objects (e.g. jewellery)
- cancer
- fever
- medical conditions, weakness
- skin abrasions or damage
- previously used fillers, Botox (6 weeks since the treatment must have passed)
- viral, bacterial and fungal skin diseases
- acute forms of rosacea and pyoderma
- active cold sores
- epilepsy
- migraines
- untreated and unregulated thyroid diseases.

Indications

- loss of firmness and elasticity of the skin
- grey, sallow complexion
- presence of wrinkles and furrows
- discolouration
- acne
- dryness and dehydration of the skin

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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The importer hereby declares that the Beautifly B-Glossy Blush device complies with Directive 2014/35 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use in certain voltage limits and Directive 2014/30 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility. The full text of the declaration of conformity is available at the following internet address: (link below)

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