BEAUTIFLY

B-Longa

PRO

series

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Laser for hair growth

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User manual

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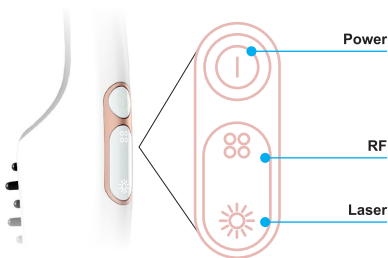
Laser comb takes advantages of the RF and laser technologies to ease the hair problem in the principle of optical hair care and high frequency hair care.

Features

- 10 beans of laser for stimulating the hair follicles.
- RF technology for stimulating scalp.
- Titanium alloy for ESD purpose, free from tying up hair.
- Simple operation with vibration + voice prompt.
- Base design for preventing scratching the product.
- Intelligent timing function to allow scientific hair care in 15 minutes.

Applicable to

- Hair loss.
- Alopecia areata.
- Androgenic alopecia.
- Hair weakened due to eg illness, narcosis, pregnancy.

Usage

1. Laser for B-Longa hair growth should be used according to the action schedule, 15 minutes, 3 times a week, preferably every other day, eg Monday, Wednesday, Friday. It is recommended to use for another 6 months.
2. The skin on the head should be freshly washed and the hair dry.
3. B-Longa should be used before applying styling products, as they can hinder the flow of laser light.
4. Press and hold the power button to start the device. Then press the RF function of the radio wave  and the Laser  to start the treatment. Place the B-Longa at the base of the hair and move from the top of the head to the nape (by combing hair). Repeat this operation for 15 minutes. The automatic timer will automatically turn off the device after the treatment.
5. After the treatment, clean the head with a cotton, soft, slightly damp cloth. Do not use strong cleaning detergents on any part of the device.

Operation instructions

1. Power on: press and hold power key  to power on, and the power key indicator will be on.
2. RF function: in power on state, press RF key  to enable the RF when the RF indicator  will be on; press again to disable the RF.
3. Laser function: in power on state, press laser key  to enable the laser when the laser indicator will be on; press  again to disable the laser function.
4. Power off: in power on state, press and hold the power key  to power off the laser comb, and the power indicator will be off.
5. Charging: connect the attached USB cable to adapter power supply to charge the laser comb. The power key will flash in red during charging, and on in red when it is fully charged. (For safety purpose, all functions of the comb do not work during charging) .

Storage and protection

Cleaning and maintenance

1. When the comb is dirty, use dry soft cloth or towel to wipe it.
2. After the treatment, clean the head with a cotton, soft, slightly damp cloth. Do not use strong cleaning detergents on any part of the device.
3. It is strictly forbidden to leak the water or other liquid into the comb.


Safekeeping

1. Keep out of the reach of children.
2. Do not put it in a place with direct sunlight, high temperature and humidity.
3. Keep it in a dry and well-ventilated place.
4. Do not disassemble, repair or modify the product without permission.

Attention

1. Do not directly stare at the red light or the red light in the mirror, which may irritate your eyes temporarily. In addition, do not use magnifier when you're using the laser comb, otherwise it may cause temporary irritation to eyes.
2. Do not put laser comb in water so as to prevent electric shock. No warranty is provided for any man-made damage to the laser comb.
3. Do not use laser combs near any heat source accessory, which may lead to failure of normal operation.
4. Laser -Avoid direct eye contact.
5. If you're taking any medicines or products that may cause you to be afraid of light, please consult the doctor first.
6. Care should be taken for women in suckling period or menstrual period.
7. Do not disassemble, modify or transform the product without permission.
8. Do not hit the product because it may cause internal damage may be caused.
9. Pregnant woman, people with sensitive skin, heart disease, abnormal blood pressure, cancers, patients with cerebrovascular diseases, acute diseases, or others receiving doctor's therapy are not allowed to use this product.
10. Disposal of the product shall be subject to the related regulations of the city on environmental protection. Otherwise, pollution may be caused to the environment.

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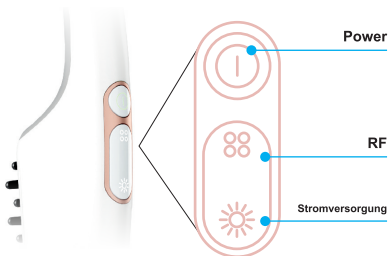Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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The importer hereby declares that the Beautifly B-Longa device complies with Directive 2014/35 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use in certain voltage limits and Directive 2014/30 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility. The full text of the declaration of conformity is available at the following internet address: (link below)

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