BEAUTIFLY

B-Lumix
Body
series

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

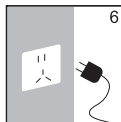
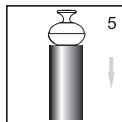
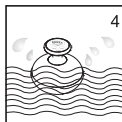
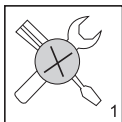
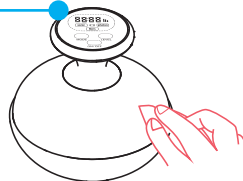
Equipment of the device

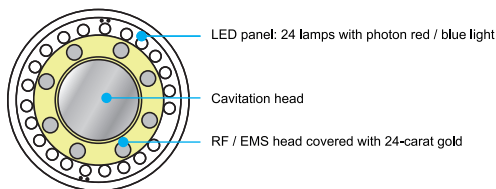Display and control panel - description

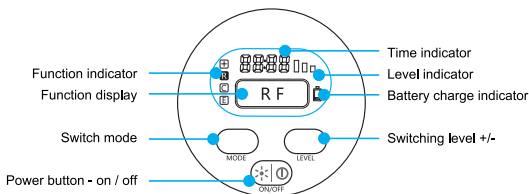Technical data B-Lumix

Input:	110-240V
Output:	5V 1A
Lithium battery:	900mAh
Dimensions (mm):	99 * 99 * 94
Cavitation frequency:	800kHz
Phototherapy lights:	red / blue
Turn-off time:	after work 5 minutes

Precautions

Before using the B-Lumix slimming and firming body device, please read the following basic safety information carefully. Observe the following precautions when using the device:

1. If your skin is very sensitive, shorten the operating time the first time you use it.
2. Avoid places where the device may fall into a sink or bathtub. If the device has been dropped in water, do not reach for it before unplugging the power cord.
3. Do not use any charger other than the original one to charge the device.
4. Do not use the device if it is damaged, wet or not working properly.

5. Do not use the device if you are pregnant or menstruating, have a pacemaker, suffer from epilepsy, heart problems or severe skin diseases.
6. Do not use the device on an area of skin that has undergone plastic surgery.
7. Do not use the device on the mouth, eyeballs or throat.
8. The device is intended for use on healthy skin. Do not use on healing or open wounds.
9. Please stop using the device if you feel uncomfortable and consult it immediately with your doctor.

B–Lumix device functions

The slimming and firming body device B–Lumix with a head covered with 24-carat gold is used for:

1. Silhouette line modeling
2. Reduction of excess body fat
3. Eliminating cellulite
4. Firming, increasing elasticity and improving the smoothness of the body
5. Skin lifting

Radio waves (RF):

Radiofrequency (RF) is a proven technology used in the medical and cosmetic industries. Under strictly controlled temperature, the magnetic pulse acts directly on the cells, generating heat energy. As a result, body rejuvenation using radio waves (RF) is the best non-surgical facelift for your skin. The use of radio waves will help you successfully model a shapely figure. You will even eliminate stubborn excess fat and cellulite. Radio wave stimulation goes to both the superficial and deep skin tissues. This will make your skin properly oxygenated, nourished and obtain proper microcirculation. As a result, along with the reduction of unnecessary fat tissue, you have the body firmness ensured all the time.

Cavitation (CAVI):

Designed for all skin types. Thanks to the strictly controlled frequency, the cavitation penetrates the skin and vibrates at a very high frequency (800KHz). Functions: mechanical, chemical and heating massage even deep skin cells, effectively shrink fat cells and increase the permeability of the cell membrane. In this way, you will improve circulation, cleanse the lymphatic system and repair dull skin tone. Thanks to cavitation, you will achieve amazing skin cleansing effects. They are visible much faster than with other treatments, and the result of a beautiful, smooth and youthful body lasts for a long time.

At the same time, cavitation will help your skin increase the absorption of personal care products. At the same time, you can be sure that the cavitation treatment is painless and completely safe.

EMS:

EMS electrostimulation stimulates the muscles to work. When using the device, the muscles contraction is identical to that during exercise. However, 50 times more energy is used than with regular exercise. Such muscle stimulation will allow you to maximally engage the muscles and increase the level of protein and elastic fibers. As a result, you will radically reduce the level of excess body fat and cellulite. EMS will restore your skin's elasticity, elasticity and guarantee your body firmness.

Red photon light:

Significantly stimulates the production of collagen proteins. Massage with the use of red photon light will allow you to visibly rejuvenate the body by improving its smoothness and elasticity. During the reduction of adipose tissue, you will radically stimulate your bloodstream and thus avoid sagging skin. Using red photon light will give your body vitality, firmness and health. This type of light is great for reducing scars, stretch marks and cellulite, as well as reducing and alleviating inflammation.


Blue Photon Light:

Especially suitable for sensitive skin. It has a positive effect on reducing inflammation, destroying bacteria and soothing irritation. Blue photon light restores the proper functioning of the skin. Its consistent use supports the reduction of cellulite and helps to balance the oil and hydration of the skin.


Charging

1. Before first use, charge the device for at least 3 hours to activate the internal lithium battery.
2. No power: When the battery runs out, the device beeps 10 times and then turns off.
3. Charging: When the device is placed on the charging cradle, the charging indicator will appear on the screen and the battery icon will gradually change over and over again - from empty to full. The backlight turns off after 30 seconds of charging.
4. Charging complete: When charging is complete, the charging indicator disappears and the battery icon shows full.

Display while charging

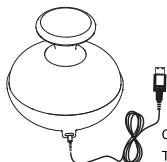Charging time: approximately 3 hours
Total working time: about 1.5 hours

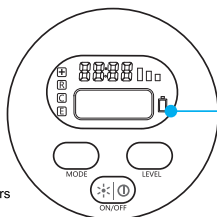Battery icon / Charging indicator


Before use

Confirmation that the B-Lumix device is ready for cavitation treatment:

Place the cavitation head flat (head facing up), select the CAVI function and drop 1-2 drops of water on the head. The water droplets will vibrate and be sprayed by cavitation.

Allergy test:

1. Clean one area of the body on the arm.
2. Apply a skin care product (gel / lotion) to cleansed skin, turn on the device and select the function you prefer. Massage this area for 2-3 minutes as instructed in this manual, then leave this area of your body without any additional treatments.
3. Wait and check for any irritation or redness within 24 hours.

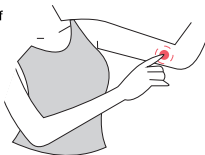User manual for combination of functions:

1. Cleanse the skin.
2. Evenly apply the skin care product (gel / lotion) to the skin.
3. Turn on the device.
4. Press MODE once to run COMBI function. Choose your preferred intensity level, RF, Cavitation and EMS will run for 10 seconds and change over and over again.
5. **Tip:** For best results, always massage your body with small circles. If necessary, add a gel / lotion during the treatment.

Radio Frequency (RF) Operation Manual

1. Cleanse the skin.
2. Apply the skin care product (gel / lotion) evenly to the skin.
3. Turn on the device.
4. Press MODE twice to start the RF function. Choose your preferred intensity level.
5. **Tip:** For best results, always massage your body with small circles. If necessary, add a gel / lotion during the treatment.

Cavitation Operation Manual (CAVI)

1. Cleanse the skin.
2. Evenly apply the skin care product (gel / lotion) to the skin.
3. Turn on the device.
4. Press MODE three times to start the CAVI function. Choose your preferred intensity level.
5. **Tip:** For best results, always massage your body with small circles. If necessary, add a gel / lotion during the treatment.

EMS operating manual

1. Cleanse the skin.
2. Evenly apply the skin care product (gel / lotion) to the skin.
3. Turn on the device.
4. Press MODE four times to start the EMS function. Choose your preferred intensity level.
5. **Tip:** For best results, always massage your body with small circles. If necessary, add a gel / lotion during the treatment.

Massage instruction

For best results, massage as follows:

Both upper arms

Massage both upper arms evenly, giving each arm the same amount of time.

Massage the upper arms from the elbows toward the armpits.

Massage the upper arms with small circles.

This form of massage relieves muscle tension, shapes the body, firms and smooths the skin.

Waist and belly

Massage your waist up and down.

Then massage the waist from the outside in.

Make a circular massage around the navel (clockwise).

Then massage the abdomen with large circles (also clockwise).

This massage of the waist and abdomen will help you quickly break down even the most stubborn adipose tissue, thanks to which you can easily get rid of its excess.

Leg massage

Massage from the ankles up.

Massage your legs in straight lines up for a while (always towards the heart).

Then massage with small circles - from the bottom of the leg up (towards the heart).

By focusing on the parts of the body particularly susceptible to cellulite, you will effectively get rid of body fat. At the same time, you can be sure that your body will not lose its firmness.

Massage of the hips and buttocks

Massage both parts of the body with movements from the bottom up.

The massage performed in this way serves to lift the buttocks and firm the skin.

Foot massage

Massage the bottom of your feet with an equal amount of time.

Massage the feet in a spiral motion.


This way you will be relieved of foot swelling. You will improve circulation, you will feel relaxed and light - you will get the effect of rested legs.

Warnings

- Treatments should not exceed 10 minutes for one body part and should be performed no more than 3 times a week. The time between each treatment should include at least one day without treatment.
- When using COMBI, RF, CAVI or EMS functions, always keep the right amount of care agent (lotion / gel) between the massage heads and the skin, especially when using the RF function.
- Always move the massage heads slowly during the treatment. Prolonged use of the device in one position may cause skin irritation.





After using the device

1. Turn off the device.
2. The skin care agent left on the massage heads, used during the massage, wipe off with a damp tissue or cloth. (Do not clean it with water and / or any chemical cleaner.)
3. Wipe and clean the device with a dry tissue or cloth. Store them in a clean and dry place.

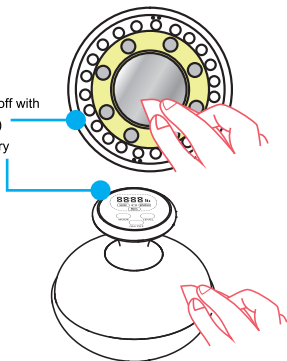
Maintenance precautions

1. Never try to disassemble the device.
2. Keep them out of the reach of children.
3. Keep the device away from fire.
4. Keep the device away from moisture. Do not submerge it in water.
5. Do not drop it or hit it.
6. Disconnect the power adapter after the device is charged.

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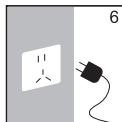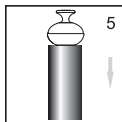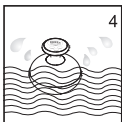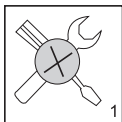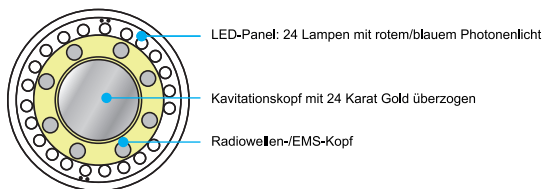Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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English

Hereby, Beautifly Sp. z o. o. declares that the device equipment Beautifly B-Lumix is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.
In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@beautifly.eu

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I represent that I have read and accepted the terms and
conditions specified in the Warranty Certificate.