BEAUTIFLY

EN

PL

B-Scrub

Blush

series

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B — Scrub
Blush
series

CAVITATION
SONOPHORESIS
PULSATING VIBRATIONS
ULTRASOUNDS
2 MODES OF OPERATION
TRANSPARENT CASE
BATTERY POWER SUPPLY

Your reliable skin scrubber. Deep hydration combined with ultrasound massage. The B-Scrub device exfoliates dead skin and cleanses the skin leaving it lighter, more flexible and healthy. Vibrations with high frequency 24 kHz soften the skin and open blocked pores. Peeling works through the use of the phenomenon of cavitation or change in the state of liquid molecules that explode under the influence of ultrasonic massage by knocking out dirt, releasing the skin from impurities and the horny layer of the epidermis.

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
1. Stainless steel shovel
2. Peeling
3. Sonophoresis
4. Switch

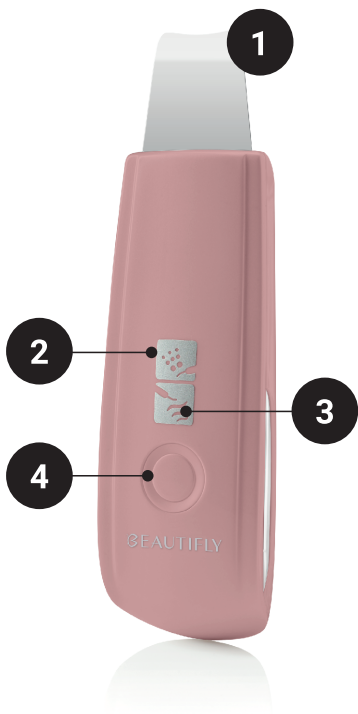Technical data

charging current:	0.3A
output current:	0.35 mA
working voltage:	DC 3.7V
charging voltage:	DC 5V
nominal voltage:	220V AC → 5V DC
work intensity:	200mA
dimensions:	4,5x18x5,5cm
vibration frequency:	24 kHz ±20%
color versions:	blush

Functions

A skin care device, uses ultrasound technology and ionophoresis, which makes it possible to exfoliate the skin, supplying nutrients to skin cells, which makes it more elastic and improves firmness.

Cavitation:

The device exfoliates dead skin and cleanses the skin leaving it lighter, flexible and healthy. Vibrations with high frequency 24 kHz soften the skin and open blocked pores. Peeling works through the use of the phenomenon of cavitation or change in the state of liquid molecules that explode under the influence of ultrasonic massage by knocking out dirt, releasing the skin from impurities and the horny layer of the epidermis.

Sonophoresis:

Sonophoresis is a procedure that uses ultrasound. Ultrasound (high frequency sound waves) increases the permeability of cell membranes. In this way, they help the nutrients contained in cosmetic preparations to penetrate into the deep layers of the skin. Ultrasounds also have a positive effect on lymph flow and cell metabolism.

Warning

- 1 Do not use this device for purposes other than those mentioned in the user manual.

Turn off the device if it is not used. Do not use the device for a

2. long time in one place. Moving the device too fast will affect its effectiveness.

3. Take away the use of the device for people with acute, chronic illness, fever, heart disease, pregnancy, facial surgery, skin rash, blood disease, during menstruation, people with pacemakers, orthodontic accommodation.
4. Do not use the device around the pupil, eyelid, injured skin, mouth and intimate parts of the body.
5. Please protect the device against physical impact.
6. Please use the device in accordance with these instructions.
Damaged device, have it repaired.
7. Please keep the device away from dust, shocks, hot or excessively cold environments.
8. This device is not waterproof, do not use it in the bathroom or in places that are overheated.
9. If the packaging is damaged, before using it, please make sure that the device is working properly.

Storage

1. Keep the device dry and clean.
2. Keep away from high temperature.
3. Keep out of the reach of children.

Cleaning

Wipe the scrub pad with a wet cloth. Do not immerse the device in water – it is not waterproof.

Before using the device

The device must be charged before the first use.

The way of charging the device

The device will turn off automatically after the "beep-beep-beep"

Open the plastic cover and connect the USB cable with the 5V 0.6A-1A power supply for charging.

The indicator  flashes during charging.

The indicator  will stop flashing after about 3 hours of charging.

Before usage

Before starting the treatment, make-up facial cleanser should be done and the make-up remover should be thoroughly cleaned. We recommend using the Beautifly B-Fresh sonic brush for this purpose. Then moisten the face with a cotton swab. During the procedure, the skin must remain moist all the time, which is why you need to do it repeat as often as necessary.

Usage

Exfoliation of the epidermis with ultrasound

1. Before starting the treatment, clean and moisturize the face.
2. Turn on the device by pressing the button once.
3. Gently slide the spatula from bottom to top, along face contours at a 45-degree angle.
4. The device will turn off automatically after 10 minutes.
5. Wash your face with water.

Negative ionization – introduction of nutrients

1. After cleansing the facial skin, apply a care product (cream, serum, ampoule, gel) on it.
2. Turn on the device by pressing the button twice.
3. Gently slide the spatula from the bottom up, along the contours of the face at a 45-degree angle. Thanks to this, you can thoroughly introduce the cosmetic.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beautifly Sp. z o. o. declares that the device equipment Beautifly B-Scrub is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address: serwis@beautifly.eu

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I represent that I have read and accepted the
terms and conditions specified in the Warranty
Certificate.