BEAUTIFLY

B-Scrub

PRO

series

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Multifunctional skin scrubber

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User manual

B-Scrub

PRO

series

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PRO

series

ultrasounds		EMS micro currents	
smart timer	IPX5 waterproof		photon light
iontophoresis	one button control		clear display
modern design		induction charging	


B—Scrub PRO is a professional device for cavitation peeling and skin firming. Cleans, cares, smoothes, brightens discolorations, exfoliates, corrects the color, but it's not all its advantages. B—Scrub PRO is also equipped with photon light, which strengthens the effects of both purification and injection of active ingredients deeply into the skin layers. B—Scrub PRO has the function of EMS microcurrents, which have the task of tightening and firming the skin, as well as raise its oval. The skin is stimulated to produce collagen and elastin, and you forget what it means to pass the time.

Technical data	
rated voltage:	3.7V
vibration frequency:	29 ± 5K Hz
EMS:	low 85V / high 100V
landing:	induction base <4H
timepiece:	5 min for each mode
size:	170mm x 52mm x 26 mm
weight:	134g
material:	ABS
waterproof:	IPX5
c colour:	white matt, rose gold

The product has a built-in lithium battery that only charges through the induction charging station. Do not place metal objects (except the B-Scrub PRO device) on the charging base to prevent burns.

Part description

1. Plug
2. Stainless steel spatula
3. LED photon light (blue / red)
4. Ultrasonic plan
5. Positive ion export plan
6. Anion import tips
7. EMS indicator (low)
8. EMS indicator (high)
9. Power indicator
10. ON / OFF / mode button
11. Induction charging base

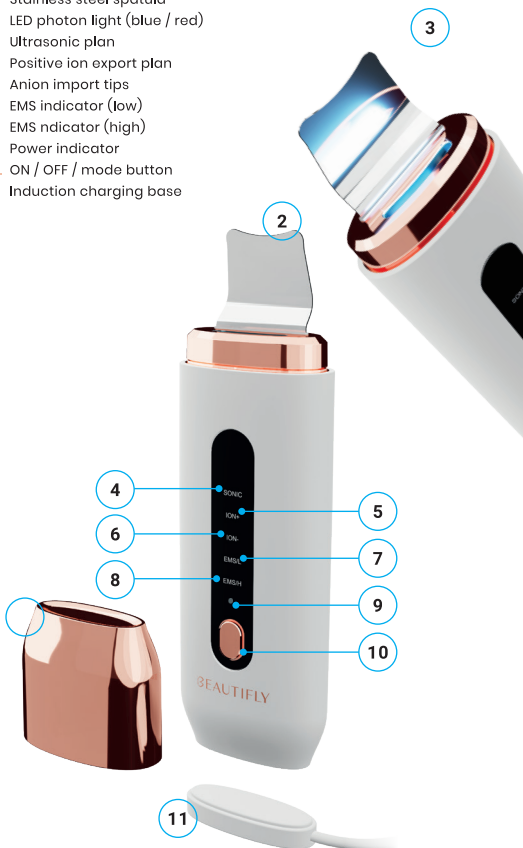Function

Multifunctional skin scrubber B- Scrub PRO for skin care. It combines ultrasound, positive ion and blue photon light technology to atomize and cleanse the dead, calloused epidermis. The combination of ultrasounds, negative ions and photon red light favors the absorption of the essence contained in used creams, ampoules or serums. Ultrasounds and EMS stimulate the production of collagen and elastin, responsible for skin tension. Apply 1-2 times a week to thoroughly cleanse the skin and improve the effectiveness of daily care.

Clean (SONIC & ION+)

After moisturizing the skin with water, high frequency ultrasonic vibration atomizes water, fully softens aging skin. Positive ion helps to export residual. Blue light prevent bacterial multiplication and balance sebum secretion.

Moisture (SONIC + ION-)

When the anion attract human micro-electricity, shovel head effectively imports the surface nutrient into basal layer. Red light with penetrability, improve the absorption of nutrient, leaving you nourishing and moisturized skin.

Anti-age (SONIC + EMS)

Microcurrents generated by ultrasound and EMS stimulate cell activity. This process stimulates the skin to produce collagen and elastin. Thanks to this, the skin becomes toned and smooth.

Description of modes

No.	Status	Indicator
1.	Clean	Power on into SONIC with ION+ mode, shovel head vibrates in high frequency, blue light on, auto-off after 5 minutes.
2.	Moisture	Power on and press 1 time on/off/mode button into SONIC + ION- mode, shovel head vibrates in high frequency, red light on, auto-off after 5 minutes.
3.	Anti-age	a. Turn on the device and press the on/off/mode button 2 times in SONIC +EMS/L mode, the shovel head vibrates at high frequency, EMS with low intensity, it automatically turns off after 5 minutes. b. Turn on the device and press the on / off button 3 times in SONIC + EMS / H mode, the shovel head vibrates at high frequency, EMS with high intensity, automatically turns off after 5 minutes
4.	Auto-off	Beeps 2 times, lights off, device off.
5.	Manual off	Press on/off/mode to skip all modes till turn it off.
6.	Low battery	Beeps 3 times, button turns red and flash.
7.	In charging	On/off/mode button turns red and keep flashing.
8.	Full charged	On/off/mode button turns green.

Usage

Exfoliation of the epidermis with ultrasounds (Ultrasounds + Ionization+)

1. Before applying the treatment, cleanse your skin thoroughly and wash your makeup. Then wet your face with a cotton swab or atomizer. For this you can use saline, tonic or water. The skin must be moist during the procedure. You cannot apply treatments on dry skin.
2. Turn on the device and start Ultrasound + Ionization + mode.
3. Gently move the spatula from bottom to top along face contours at an angle of 45 degrees. The spatula should be accessible in such a way that the photon light reaches the skin.
4. The device will turn off automatically after 5 minutes.
5. Wipe your face with tonic.

Detailed correct instructions for cavitation peeling are available on our YouTube channel:

<https://www.youtube.com/watch?v=xjL3O1BQ6A8>

Introduction of nutrients using ultrasounds (Ultrasounds + Ionization-)

1. After cleansing your face, apply the product to it care (cream, serum, ampoule, gel).
2. Turn the device on and press the On / Off switch once to activate the Ultrasound + Ionization- mode.
3. Gently move the spatula from bottom to top along the contours face. The spatula should be arranged in the opposite way than in the first step. In case of doubt about how to properly hold the device, please visit our channel on the YouTube platform.
4. The device will turn off automatically after 5 minutes.

EMS (Ultrasounds + EMS currents).

1. After facial skin care, turn on the device and press On / off switch twice to enter Ultrasound mode + EMS / weak (two work levels to choose from).
2. To improve skin elasticity, move with a spatula, by making movements from the bottom up along the contours of the face under 45 degree angle.
3. The device will turn off automatically after 5 minutes.
4. After finishing using the machine, wipe it with a damp cotton cloth.

Cautions

1. Person with electronic human organ are forbidden to use it.
2. Person with skin surgery within 1 month or skin sunburn is forbidden to use.
3. Forbidden to use by person with inflammation of the skin, scar, damaged and electromagnetic sensitive is forbidden to use.
4. Take off watch, glasses and other metal objects before operation.
5. Avoid shovel head close to eyeball.
6. Don't operate in the same area for long time.
7. Do not use the device for pregnant women.

Cleaning

Rinse the shovel head and wipe it with a dry soft cloth after using.
Do not immerse the device into water for long time.

Storage

1. Keep it dry and clean.
2. Keep it away from high temperature.
3. Keep it out of reach of children.

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The crossed-out wheeled bin symbol indicates that in the European Union after the end of the product's useful life it must be disposed of in a separate, dedicated collection point. Do not dispose of these products together with unsorted general household waste

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Beauty Factor Sp. z o.o.
Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beauty Factor Sp. z o. o. declares that the device equipment Beautifly B-Scrub PRO is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beauty Factor Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@rgmedia.pl

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